


Small in
Size
Big in **Performance**

Technical data

| Measurement range | |
|--|--|
| SpO₂ | 45 to 100 % |
| Pulse frequency | 20 to 300 BPM |
| Accuracy | |
| Saturation | +/- 2 % (70 to 100 %) |
| Pulse frequency | +/- 1 BPM to 100 BPM +/- 1% >100 BPM |
| Batteries | 2 x AA A ALKALINE |
| Capacity | >30 hrs |
| Operating temperature | 0° to + 45°C |
| Storage temperature | -20° to 70°C |
| Humidity | 0 to 90 % (non condensing) |
| Measurement process | Split pulse wave with Fuzzy Logic Control |
| Display | LCD with white backlighting and red, orange and green LED's. Indications of % SpO ₂ , pulse rate, signal quality, bar graph, battery status und sensor status |
| Pulse tone | Tone pitch variable (depending upon the measured O ₂ saturation) |
| Signal | Acoustic and visual signals (pulse frequency, sensor status, system errors and battery condition) Signal: on/off function with automatic deactivating after 3 min. |
| Dimensions | 32 x 136 x 24 mm |
| Weight | ca. 72 g with batteries |
| Classification | EN 60601-1 |
| Protection system against electric shock | Protection class II |
| Degree of protection against electric shock | Type BF |

Distributed by:



ENVITEC-Denmark ApS • Sankt Knuds Vej 37 • 1903 Frederiksberg C, Denmark
Tel. +45 43 96 43 94 • Fax. +45 43 96 43 98 • e-mail: info@envitec-dk.com • www.envitec-dk.com

The manufacturers listed are the holders of their respective names and/or trademarks, and are not to be taken as an endorsement or affiliation with Envitec- Denmark ApS.
Caution: Federal law (USA) restricts this device to sale by or on the order of a licensed physician.



O₂

OxiPen® – pulse oximeter

The most versatile
"patient to patient" pulse oximeter
for spot checks and short term monitoring.

Medical technology
Monitors



OxiPen® and SoftTip® medium



SoftTip® large



FingerClip large




SoftTip® medium



Multiseite Y- Sensor



Multiseite WRAP Sensor

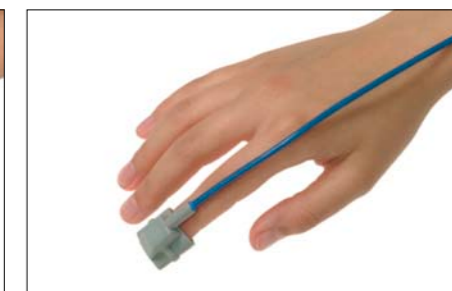

SoftTip® small

Small in size – big in performance

Proven & trusted technology from EnviteC - the world leader in SpO₂ sensor technology.

OxiPen® - culmination of years of clinical experience, brings you the latest advantages for "patient to patient" SpO₂ monitoring.

The new generation electronics and extremely handy pocket- size oximeter design provide effective solution for spot check and short term monitoring.

Whenever a quick and reliable measurement of oxygen saturation is needed, **OxiPen®** delivers fast, precise and reliable support.

While small in size **OxiPen®** is big on performance, providing SpO₂, pulse rate, signal quality and pulse strength for all patient groups.

Superior quality

OxiPen® detects permanently the quality of the pulse signal for reliable measurement and presents continuously audible and visible signal tone for the pulse rate. The LED – pulse quality indicator provides easy-to see fast assessment of patient readings.

Additionally the separate displayed **pulse-strength bar** indicates relative perfusion at the measurement site.

Thus the integrated Artifact Leveling System –**"AL"**, **OxiPen®** secures the optimal measurement results without compromising accuracy, dynamics and artifact reliability in applications.

Patient to patient monitoring

Versatility in application

The pulse oximeter **OxiPen®** offers you proven performance and ultimate versatility in application in EMS, hospital or clinical environments as well as in home use.

From pre and post operative monitoring, pulmonology, movement therapy, home care and transportation and everywhere in -between- the **OxiPen®** with its simple functionality with its small and ergonomic styling makes it a reliable partner for you.

The small and lightweight design makes it possible to carry **OxiPen®** anywhere in your pocket.

Scalable & flexible SpO₂ sensor solutions

The special developed SpO₂ sensors for **OxiPen®** accommodate a wide range of extremities from **pediatric to adult** patients and meet your application needs while keeping high **hygienic standards**.

The **OxiPen®** is available with the following SpO₂ sensors:

- **large finger SoftTip®**
- **medium finger SoftTip®**
- **small finger SoftTip®**

To secure high hygienic standards and eliminate cross contamination risks during "patient to patient" monitoring, the fully encapsulated **SoftTip®** SpO₂ finger sensors make routine cleaning and disinfection easy and safe.

The **SoftTip®** SpO₂ finger sensors are compatible with most institutional cleaning solutions, including soap, alcohol and commercial disinfectants. You can choose the protection level you require.

EnviteC's quality SpO₂ sensors offer proven performance and well known durability in usage.

Easy-to-use


User-friendly two button operating interface coupled with large illuminated display with clearly presented information secure comfortable support for the operating regardless of application needs.

OxiPen® provides beside sensor connection and low-battery indicator also automatic deactivating function after 3 minutes when not in use. It allows to provide 30 hours of continuous operation or approximately 3000 Spot checks by average duration of 1 minute.

The proven performance combined with simple operation makes **OxiPen®** your reliable guide in "patient to patient" SpO₂ monitoring.



OxiPen®

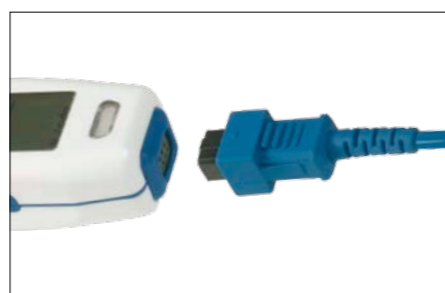

SoftTip® large



Easy cleaning and disinfection