


ENVITEC

by Honeywell

Monitors

Versatile and Precise

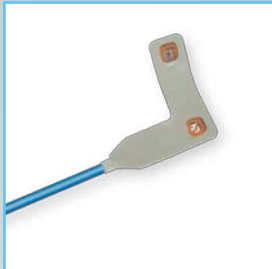

SpO₂

PulsQuant P Digital pulse oximeter

- Easy-to-use
- Compact and attractive design
- Full alarm system

Medical technology
Monitors

PulsiQuant P digital pulse oximeter

Description

PulsiQuant P enables the measurement of the arterial blood's oxygen saturation. The graphical LC-display indicates the plethysmogram curve, the arterial blood's functional oxygen saturation (SpO₂) and the pulse rate.

- lightweight, compact and attractive design
- easy-to-use menu
- SpO₂, pulse rate and pulse strength measurement
- plethysmogram, alternatively numerical screen
- integrated rechargeable battery
- designed for infants and adults
- full alarm system
- trend display for pulse and SpO₂
- memory for SpO₂ and Pulse rate for up to 12 hours
- serial Interface for printing and PC communication

Technical data

Measuring range	
SpO ₂	0 -99%
Pulse	30-249 1/min
Accuracy	
SpO ₂	+/- 2 % (70 to 99%) not specified for lower values
Pulse	+/- 1% or +/- 1 digit
Display	Oxygen saturation (actualization of the measured values with each pulse), pulse, and also plethysmogram curve (indication not proportional to the size of the pulse wave)

Patient module protected against defibrillation.

